

Dri-Stor MAINTENANCE PLANS

LogiQuip is committed to customer satisfaction throughout the life of your products. With this in mind, we are pleased to offer Maintenance Plans to keep your new LogiQuip Dri-Stor Drying Cabinet running like new. These hassle-free plans take the maintenance of your Dri-Stor Drying Cabinets off of your to-do list.

STANDARD

WARRANTY

LogiQuip Dri-Stor Cabinets are covered by a **5-Year Limited Warranty** from date of receipt against defect in materials or workmanship under normal use and service, excluding electronic components. During this period, LogiQuip will, in its sole discretion, repair defects in materials or workmanship that existed when the product was received or replace the product. If the product is replaced, the new product will be warranted for the remaining term of the original product warranty period.

1-Year Warranty on electronic components and sub-assemblies.

2-YEAR

MAINTENANCE + STANDARD WARRANTY

5-YEAR LIMITED WARRANTY + 24 MONTHS OF MAINTENANCE

- Maintenance Visits Every 6 Months
- Replacement and Consumable Parts (1) Main Cabinet Filter (40000MFF)
- Inspect Tubing
- Check Operation of Fans and Pumps
- Check Doors and Locks
- Check Levelness of Cabinets

4-YEAR

MAINTENANCE + STANDARD WARRANTY

5-YEAR LIMITED WARRANTY + 48 MONTHS OF MAINTENANCE

- Maintenance Visits Every 6 Months
- Replacement and Consumable Parts (2) Main Cabinet Filter (40000MFF)
- Inspect Tubing
- Check Operation of Fans and Pumps
- Check Doors and Locks
- Check Levelness of Cabinets

MAINTENANCE SCHEDULE

TYPE OF SERVICE	MODEL NUMBER	6 MONTHS	12 MONTHS
Cabinet Pre-Filter	40000EF	Replace	Replace
Channel Purge HEPA Filter Cannister (if equipped)	40000CPHF	Replace	Replace
HEPA Pre- and Post-Filters (including O-Rings)	40000PPFK	Replace	Replace
Scope Cabinet HEPA Filter*	40000MFF	Replace Every 2 years	
Tubing [†]	Model # Varies	Inspect	Inspect

*Can be ordered separately

[†]Tubing replaced is an additional cost



MODEL # 40000MFF

Dri-Stor MAINTENANCE KITS

Keep your Dri-Stor Drying Cabinet running like new with our Filter Maintenance Kits. Choose from three (3) Kits (below) for a one-year supply of Filters and O-Rings which need to be replaced every six months.

Kit 1

1-YEAR SUPPLY

Replace Filters and O-Rings every six months


FOR CHANNEL PURGE CABINETS WITH 10 SCOPES/PROBES OR LESS

Model 40000CPKIT-1

INCLUDES:

40000CPHF - 2 Channel Purge Hepa Filter Canisters

40000HEPAKIT - 2 HEPA Cabinet Pre-Filters, 2 HEPA Pre-Filters, 2 HEPA Post-Filters, 2 O-Rings

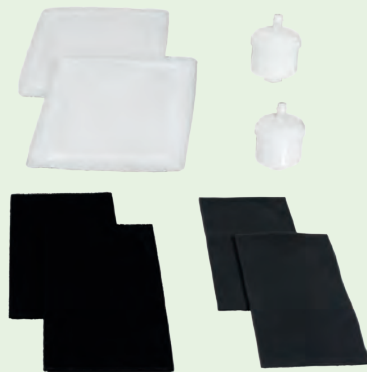

Kit 2

1-YEAR SUPPLY

Replace Filters and O-Rings every six months


FOR CHANNEL PURGE CABINETS WITH 11-20 SCOPES

Model 40000CPKIT-2

INCLUDES:

40000CPHF - 4 Channel Purge Hepa Filter Canisters

40000HEPAKIT - 2 HEPA Cabinet Pre-Filters, 2 HEPA Pre-Filters, 2 HEPA Post-Filters, 2 O-Rings

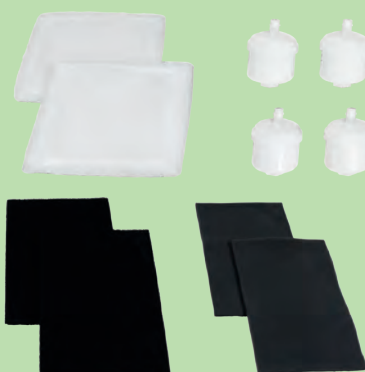

HEPA Kit

1-YEAR SUPPLY

Replace Filters and O-Rings every six months

FOR HEPA FILTER MODELS WITHOUT CHANNEL PURGE

Model 40000HEPAKIT

INCLUDES:

2 HEPA Cabinet Pre-Filters

2 HEPA Pre-Filters

2 HEPA Post-Filters

2 O-Rings

