



Case

CARTS



TABLE OF CONTENTS

TABLE OF CONTENTS

LOGIQUIP® CASE CARTS

FEATURES & BENEFITS 3-4

STAINLESS STEEL BOLTED CASE CARTS 5

SHELVING OPTIONS FOR BOLTED CASE CARTS 5

OVER-THE-ROAD CASE CARTS, ACCESSORIES & CASTERS 6

CASE CARTS ACCESSORIES & CASTERS 7



LOGIQUIP® CASE CARTS

Designed specifically to protect sterile supplies as they move from central processing to surgical suites, LogiQuip **Case Carts** provide a safe, durable, and efficient means of transportation. All **Stainless Steel Bolted** and **Over-the-Road Case Carts** have a full wrap-around **Bumper**, **6" Casters** for easy maneuverability, and accept height adjustable **Shelving** (**Shelves** not included). LogiQuip also offers **Open Case Carts** (custom option) for a lighter option that is the ideal solution for a dedicated clean or soiled delivery system from the OR to reprocessing.

STAINLESS STEEL CONSTRUCTION

An easy-to-clean, corrosion-resistant Stainless Steel body on all **Closed** and **Over-the-Road Case Carts** provides strength and durability in the harshest environment.

6-INCH CASTERS

These durable **Casters** can withstand high temperatures and high pressure when going through an automated washer. Sealed ball bearings allow these **Casters** to quietly move through your facility.

ACCESSORIES & COMPONENTS

With a variety of bolt-on accessories and components including **Card Holders**, **Clean/Dirty Wheel Indicator**, **Handles**, and **Bumpers**, you can customize these **Case Carts** to help streamline processes.

FULLY SERVICEABLE CONSTRUCTION

Case Carts are abused daily, which can lead to damage. With **Bolted** (vs. welded) **Construction**, LogiQuip **Case Carts** are easily serviced to reduce downtime and extend the life of your **Case Carts**.



OVER THE ROAD CASE CART



BOLTED CASE CART



OPEN CASE CART
(Available as Custom)

LIMITED 5-YEAR WARRANTY

LogiQuip products are covered by a **Limited 5-Year Warranty** from date of receipt against defect in materials or workmanship under normal use and service, excluding **Electronic Components**, **Cart Covers**, **Custom Products**, **LogiCell Tambour Doors**, **Label Holders**, and **Casters**. During this period, LogiQuip will, in its sole discretion, repair defects in materials or workmanship that existed when the product was received or replace the product. Visit LogiQuip.com/Warranty for more details.

LOGIQUIP® CASE CARTS



LOGIQUIP BOLTED STAINLESS STEEL CASE CARTS

BOLTED STAINLESS STEEL CASE CARTS

LogiQuip **Bolted Case Carts** are easily serviceable, reducing downtime and extending the life of your **Case Carts**. Large enough to accommodate oversized instrument sets, these **Case Carts** come standard with **Ergonomic Handles**, **270° Door Swings**, **6" Casters**, and **Wrap Around Bumper**. The **Stainless Steel Cart Body** is easy-to-clean, has no water retention points, and comes with a no drip "tilted" pan base. **Shelving** options are sold separately.

BOLTED STAINLESS STEEL CASE CARTS		
MODEL	EXTERIOR D x W x H - IN.	INTERIOR D x W x H - IN.
LQCC3630	27 x 36 x 30	25 x 30.5 x 22
LQCC3640	27 x 36 x 40	25 x 30.5 x 32
LQCC3655	27 x 36 x 55	25 x 30.5 x 47
LQCC4640	27 x 46 x 40	25 x 40.25 x 32
LQCC4655	27 x 46 x 55	25 x 40.25 x 47
LQCC5240	28 x 52 x 40	25 x 46.5 x 32
LQCC5255	28 x 52 x 55	25 x 46.5 x 47




MODEL # LQCC3640

SHELF OPTIONS FOR BOLTED CASE CARTS

There are three (3) Roll-Out **Shelving** options available for **Bolted Stainless Steel Case Carts**; **Perforated Stainless Steel Shelves**, **Solid Stainless Steel Shelves**, and **Wire Case Shelves**. All **Shelves** are sold separately and come with **Adjustable Rails**. Maximum weight capacity is 200 pounds for **Stainless Steel Shelves** and 150 pounds for the **Wire Shelves**.


BOLTED CASE CART SHELF OPTIONS		
MODEL	DESCRIPTION	D x W - IN.
LQCCS36PR	Roll-Out Perforated Stainless Steel Case Cart Shelf for LQCC3630, LQCC3640, and LQCC3655	23 x 26
LQCCS36SR	Roll-Out Solid Stainless Steel Case Cart Shelf for LQCC3630, LQCC3640, and LQCC3655	23 x 26
LQCCS36WR	Roll-Out Wire Case Cart Shelf for LQCC3630, LQCC3640, and LQCC3655	23 x 26
LQCCS46PR	Roll-Out Perforated Stainless Steel Case Cart Shelf for LQCC4640 and LQCC4655	23 x 36
LQCCS46SR	Roll-Out Solid Stainless Steel Case Cart Shelf for LQCC4640 and LQCC4655	23 x 36
LQCCS46WR	Roll-Out Wire Case Cart Shelf for LQCC4640 and LQCC4655	23 x 36
LQCCS52PR	Roll-Out Perforated Stainless Steel Case Cart Shelf for LQCC5240 and LQCC5255	23 x 42
LQCCS52SR	Roll-Out Solid Stainless Steel Case Cart Shelf for LQCC5240 and LQCC5255	23 x 42
LQCCS52WR	Roll-Out Wire Case Cart Shelf for LQCC5240 and LQCC5255	23 x 42

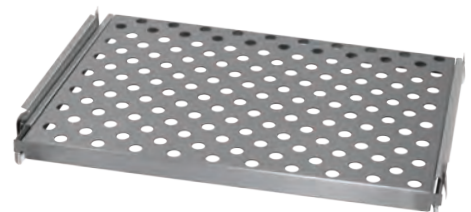

Roll-Out Perforated Stainless Steel Shelf



Roll-Out Solid Stainless Steel Shelf

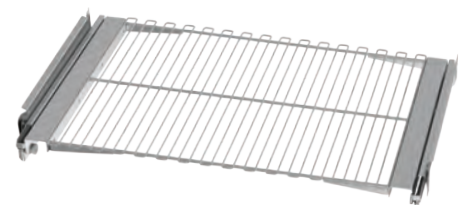

Roll-Out Wire Shelf

OVER-THE-ROAD CASE CARTS

OVER-THE-ROAD CASE CARTS

Building on our **Standard Case Carts**, LogiQuip's **Over-The-Road Case Carts** provide a rugged, durable and safe way to transport sterile supplies when Central Processing is off-site. All welded, stainless construction, heavy-duty 6" x 2" washable casters and optional accessories including a **Tow Hitch**, help to create a clean and safe environment.

- Exterior Dimension: 30" D x 43.75" W x 57" H (8.28 sq. ft.)
- Shelf Size: 26" D x 39" W
- Capacity: Six (6) full-size containers
- Container Capacity: Available for 12, 18, and 24 containers
- Unique container retention system
- Storage Space: Top 5"H; Bottom 14.5" H
- Casters: Two (2) directional locking; Two (2) swivel with brake
- Bumper System: Full wrap around on top and bottom of cart
- Door System: Overlapping, self-sealing with 95° cart seal and double locking key
- Security Tab: Breakaway



MODEL #
OTRCC3657SS



MODEL # OTR0100-01-KT Easy slide
CLEAN/DIRTY Indicator eliminates
guesswork and errors



MODEL#
CW31-76RXSR
Rigid Caster

MODEL#
CW31-76BXS
Swivel Caster w/Brake

OVER-THE-ROAD CASE CART CASTERS/ACCESSORIES

MODEL	DESCRIPTION
OTR0100-01-KT	Clean/Dirty Indicator
CW31-76RXSR	Cart Caster – Rigid
CW31-76BXS	Cart Caster – Swivel with Brake
CC-FLAG3	Flag System
OTRCC-TOW	Over-the-Road Case Cart Tow Hitch



MODEL# OTRCC-TOW
Tow Hitch System



MODEL# CC-FLAG3
Flag System

CASE CARTS ACCESSORIES & CASTERS

CLOSED CASE CART ACCESSORIES

Accessories for LogiQuip's Closed Case Carts are shown at right and below:

CLOSED CASE CART ACCESSORIES	
MODEL	DESCRIPTION
BO-0204	Pick List Card Holder – 8.5" x 11" Stainless
CC-TSE	Top Extension Drop Shelf – 16" x 24"
CC0230-01-BIO	CLEAN/DIRTY Indicator Wheel
CC-TOW	Tow Hitch System
CC-FLAG3	Flag System (shown on page 6)



MODEL# CC0230-01-BIO
CLEAN/DIRTY Indicator Wheel



MODEL# CC-TSE
Extension Drop Shelf



MODEL# BO-0204
Card Holder



CC-TOW
Tow Hitch System

CLOSED CASE CART CASTERS

Our **Standard Case Cart Casters** are 6" x 1.5" stainless steel, all swivel (one is directional locking and two are braking), and have the following:

- Sealed bearings - upper and lower
- Grease fitting - bearings and axle
- High temp wheels - poly
- Cart washable
- Optional sizes available upon request



MODEL#
CCCAST6



MODEL#
CCCAST6B



MODEL#
CCCAST6DL

CLOSED CASE CART CASTERS	
MODEL	DESCRIPTION
CCCAST6	Cart Caster – Swivel
CCCAST6B	Cart Caster – Swivel with Brake
CCCAST6DL	Cart Caster – Directional Lock



www.logiquip.com

(800) 665-3760 | customersupport@logiquip.net

9.2022