

NEW 

DECONTAMINATION

SINKS – FIXED & HEIGHT ADJUSTABLE

LogiQuip's **NEWLY REDESIGNED Decontamination Sinks** are more ergonomic than ever! NOW offering four (4) **Height Adjustable Options** that adjust with the press of a button, our **Decontamination Sinks** now come standard with 20" long **Basins**, a 12" **Backsplash** that allows for a **Universal Workflow**, **Marine Edge Countertops**, and **Basins** that are positioned closer to the user for less reaching. Add the optional **Pegboard Backframe** and **Optional Accessories**, including **Overhead LED Lighting**, to create an ergonomic and efficient workspace.

HEIGHT ADJUSTABLE

Easily raise and lower with the push of a button. Our **Electronic Height Adjustable Decontamination Sink** adjusts from 35" to 42" high, providing an ergonomic workspace for any user. **Fixed Height** options are also available.

300 SERIES STAINLESS STEEL

LogiQuip **Decontamination Sinks** are **300 Series Stainless Steel** making them durable, easy-to-clean* and corrosion-resistant. All welded construction creates virtually seamless connection points that minimize cracks and crevices where germs can hide.

UNIVERSAL WORKFLOW


Faucets and **Water** and **Air Guns** are mounted to the 12" **Backsplash** with openings over each **Basin**. This allows you to connect and move the **Faucets** and **Guns** to create the workflow you need. **Hole Caps** are provided to cover openings that are not used.

24-INCH BASINS

LogiQuip **Decontamination Sinks** come standard with 24" wide and 20" long **Basins** to give more space for cleaning instruments. Available in **Single**, **Double**, and **Triple Basin** configuration, each **Basin** has **Etched Gallon Marks** to help ensure precise chemical dilution to safe and effective cleaning. 30" wide **Basins** are also available as an option.

OPTIONAL ACCESSORIES

Several additional options including **Pegboard Backframe** and **Accessories Options**, **Overhead LED Lighting**, **Eye Wash**, **Dosing Soap Pumps**, and **Foot Operated Drains** allow you to customize our **Decontamination Sinks** to your specific needs.

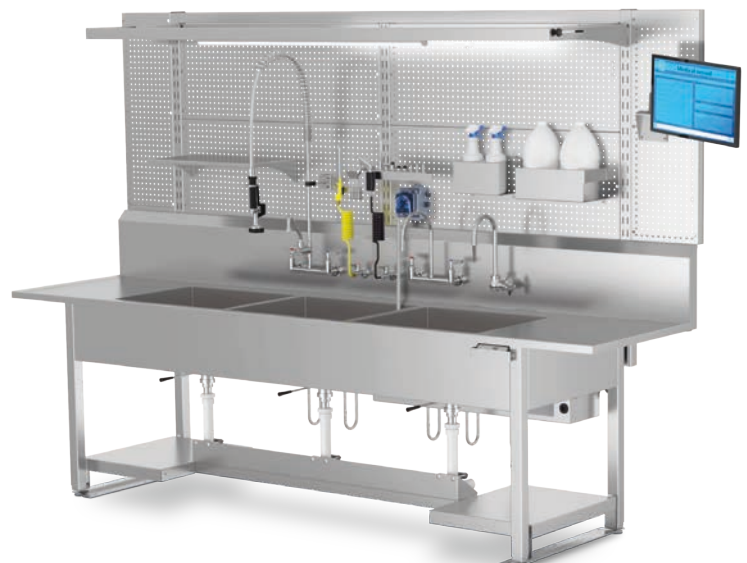

*See the LogiQuip Cleaning Guide for cleaning instructions.

LOGIQUIP® DECONTAMINATION SINKS

FIXED HEIGHT DECONTAMINATION SINKS – FEATURES & MODELS				
FEATURES	LQDSF307237 Single Basin Sink 30" x 72" x 37"	LQDSF309637D Double Basin Sink 30" x 96" x 37"	LQDSF309637T Triple Basin Sink* 30" x 96" x 37"	LQDSF3012037 Triple Basin Sink 30" x 120" x 37"
Mixing Valve Faucet(s)	1	2	2	2
Pre-Rinse Faucet(s)	1	2	2	2
Water Gun Kit (8 Tips)	✓	✓	✓	✓
Air Gun Kit (5 Tips)	✓	✓	✓	✓
Basin Size 20" x 24" x 10"	1	2	3	3
Etched Gallon Marks	✓	✓	✓	✓
12-Inch Backsplash	✓	✓	✓	✓
Wrap Around Apron	✓	✓	✓	✓
Open Pipe Base	✓	✓	✓	✓
Lever Handle Drains	1	2	3	3
Adjustable Feet	✓	✓	✓	✓
Seismic Bracing	✓	✓	✓	✓

ADJUSTABLE HEIGHT DECONTAMINATION SINKS – FEATURES & MODELS				
FEATURES	LQDSA307235 Single Basin Sink 30" x 72" x 35-42"	LQDSA309635D Double Basin Sink 30" x 96" x 35-42"	LQDSA309635T Triple Basin Sink* 30" x 96" x 35-42"	LQDSA3012035 Triple Basin Sink 30" x 120" x 35-42"
Mixing Valve Faucet(s)	1	2	2	2
Pre-Rinse Faucet(s)	1	2	2	2
Water Gun Kit (8 Tips)	✓	✓	✓	✓
Air Gun Kit (5 Tips)	✓	✓	✓	✓
Deionized (DI) Water Faucet	✓	✓	✓	✓
Basin Size 20" x 24" x 10"	1	2	3	3
Etched Gallon Marks	✓	✓	✓	✓
12-Inch Backsplash	✓	✓	✓	✓
Wrap Around Apron	✓	✓	✓	✓
Open Pipe Base	✓	✓	✓	✓
Lever Handle Drains	1	2	3	3
Adjustable Feet	✓	✓	✓	✓
Seismic Bracing	✓	✓	✓	✓

*96" Wide Triple Basin Sinks have less counterspace than 120" Wide Triple Basin Sinks.